

PRO-NOX™ Options




(Minimum) peak flow rate	120 LPM
Inhalation resistance (w/6' circuit)	0 to -7.5 cm, H2O @ 60 LPM
Initiation pressure	-2.5 cm H2O
Input Pressure	50 to 70 PSI (3.5 to 5 Bar)
Operating Temperature	41° to 104° F
Storage Temperature	-40° to 140° F
Pro-Nox™ Mixer Dimension (L x W x H)	9" x 6.5" x 3.5"
Portable E-Cylinder Mobile Cart:	47.5" w/Mixer on Stand
Overall Height Stand	44" high ~ stand only
Stand width at Base	23" wide
Weight	26 lbs. complete hospital system
Oxygen Inlet	DISS 1/2" OD Green Hose
Nitrous Inlet	DISS 3/8" OD Blue Hose
Pro-Nox Active Scavenger Vacuum Connection	1/4" vacuum hose on 1/4" hose barb
Scavenger Hose	20' long
Vacuum flow (Dependent on vacuum source)	35–50 LPM Maximum (Ball Lever in full on position)



Pro-Nox
E-Cylinder
Mobile Cart



Pro-Nox Wall Mount

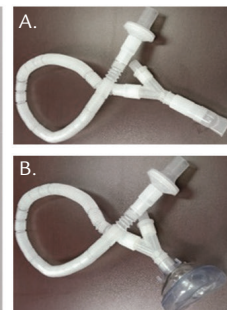

A. Patient Circuit
with Mouthpiece
B. Patient Circuit
with Mask



Pro-Vac Portable Scavenging System

PRO-NOX™

The Premier Name in Nitrous Oxide Systems



Patients using Pro-Nox in procedure rooms

Nitrous Oxide and Oxygen Analgesia System Great for any painful office procedure!

- Self-Administered by Patient
- Quick Onset
- Fixed 50% N2O & 50% O2
- Pain Management
- Short Duration of Effect



125 Fleming Drive, Cambridge ON N1T 2B8
www.clarionmedical.com



125 Fleming Drive, Cambridge ON N1T 2B8
www.clarionmedical.com

PRO-NOX™

The Newest Innovative Technology

Our team has worked to design a system that addresses all the needs of your office staff AND the patient! **Super Simple and Easy to Use**, there is no need to read cumbersome regulators or understand PSI conversion charts to guess when you need to change tanks. With the PRO-NOX™ Nitrous Oxide Delivery System, you have audible and visual alarms that warn you when either gas is running low. We offer superior infection control with our internal on-demand valve. No attached parts that need to be wiped down that have a chance to break down or get lost. All important working parts are located inside the small, state of the art mixer.

The CAREstream team is here to support you through continuous learning on our many webinars, training sessions and more! In addition to the outstanding support personnel, we can provide you with **Articles, Studies, References** and so much more. We are here to help you through the process of adding Nitrous to your pain management toolbox.



Pro-Nox Nitrous Oxide System

Benefits of PRO-NOX™ 50% O2/50% N2O Analgesia

- Proven effective & safe—used globally for decades.
- Patient administered—empowers the patient to respond to their pain & anxiety when they need it.
- Quick Onset—takes effect in seconds—immediate relief!
- Short duration of effect—clears from body within minutes.
- Additional pain tool in your toolkit that is easy to administer.
- Can be used for the duration of any procedure where the patient has pain and anxiety.
- Speeds patient recovery allowing them to drive to and from their procedure!
- Non invasive—easy to administer and set-up.
- Improves patient satisfaction scores!
- Low cost to use and minimal disposable cost.

The PRO-NOX™ Aesthetic System

comes complete with the Pro-Nox mixer, stand, E-tank holder, basket, regulators for N2O & O2, along with quick guides, five easy steps, instructional video and support for your staff.

Safety Features of PRO-NOX™

- Fixed 50/50 N2O and O2 delivery
- Internal On-Demand Valve
- Superior Infection Control
- Single-Patient use Circuit with oneway valve
- Audible and Visual low gas alarms
- Safety shut off when gas pressure too low to deliver 50/50 mix
- Patient Self-Administered
- Scavenger system available
- Ready Made Kits—Quick & Easy!

A Safe Mixture of 50% Nitrous Oxide and 50% Oxygen for pain & anxiety

The 50/50 mix of Nitrous Oxide and Oxygen has been used around the world in labor and delivery for decades. Come join the growing number of private practices and clinics around the world that offer the PRO-NOX™ Nitrous Oxide Delivery System to help ease pain and anxiety. **Patients love the Pro-Nox because they control when to use the system**, on a fully portable cart, allowing them to move between treatment rooms, while still getting the relief they need. The 50% N2O/50% O2 is an effective analgesic that provides quick pain relief and allows for quick recovery all while giving the patient a sense of control. Studies have shown that it is non-addictive for the patient for procedures, while delaying or decreasing the need for other pain medicine. From Aesthetics, Dermatology, Plastic Surgery, Vascular, OBGYN, Orthopedics, Pain Management, Regenerative Medicine, Urology, Birth Centers, Hospitals, Surgery Centers and More ~ **Pro-Nox is a great alternative for pain management.**

Patients LOVE the PRO-NOX™ System and are happy to share their experiences:

“ I couldn’t believe how easy it was to use and how quickly it took effect. It was enough relief to get me through my liposuction and cell therapy procedure.”

- Patient from Florida

Doctors & their staff LOVE the PRO-NOX™ System and are happy to share their experiences:

“ The Pro-Nox allows me to reduce the amount of other pain medicine required and enables the patient to drive to and from procedures that they would otherwise require a driver. Helps differentiate my practice and my patients love it ”

- Dr. Jason Emer ~ Beverly Hills, CA

“The Pro-Nox is the best thing I have ever bought! My patients love it and it makes the procedures go so much smoother.”

- Dr. Warren B. Seiler III ~ Homewood, AL
Executive Director ~ American Board of Laser Surgery

- | | | | |
|------------------------|-------------------------------------|-------------------------------|-----------------------------------|
| - Blood Draws | - EMS | - Labiaplasty | - Sclerotherapy |
| - Breast Augmentation | - Endometrial Biopsy | - Labor & Delivery | - Stem Cell Therapy |
| - Burn Debridement | - Endovenous Laser Treatment (EVLT) | - Laser Treatments | - Tattoo Removal |
| - Butt Lifts | - Face Lifts | - LEEP Procedure | - Thermage |
| - Circumcision | - Facial/Body Injections | - Liposuction | - Tummy Tucks |
| - Colposcopy | - Fat Transfer | - Micro-Needling/PRP | - Ultherapy |
| - Dilation & Curettage | - Hair Restoration | - Orthopedic joint injections | - Vasectomy |
| - ED Treatment | - IUD Insertion | - Permanent Make-up | - Any Procedure with Pain! |

